


## TRB-624

Column: **TRB-624**, P/N TR-601462  
 Dimensions: 60m x 0.25mm x 1.4µm  
 Injection: 1µL solvents mixture, split 1:100 (20-600 ng/comp.), 260°C  
 Carrier gas: H<sub>2</sub>, constant pressure 25 psi (172.3 KPa).  
 Oven temperature: 50°C(5min) @ 6°C/min to 220°C  
 Detector: FID, 280°C

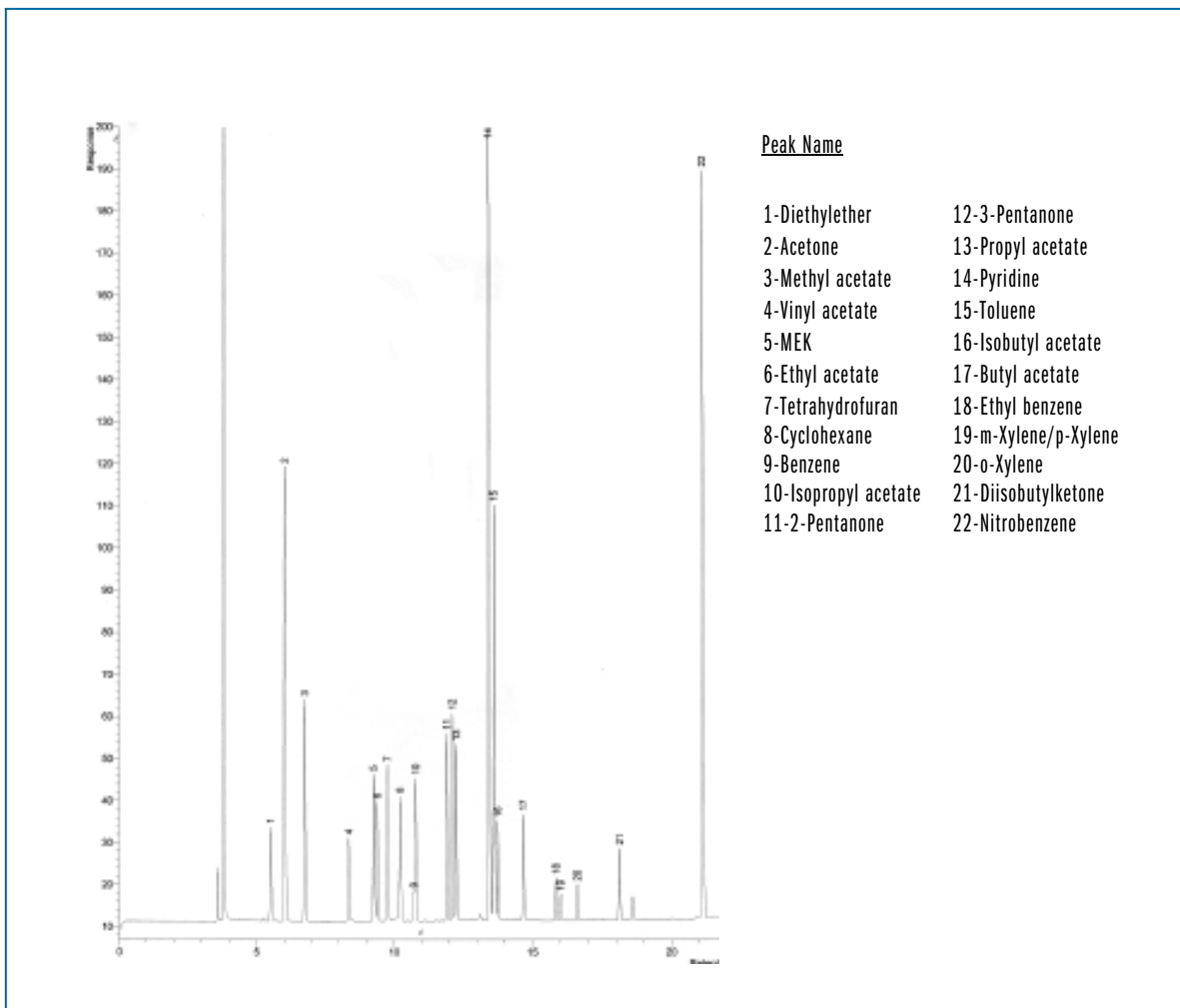

TKG 1119

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es