

## MIXTURE OF SOLVENTS AND ISOMERS OF N-HEPTANE

Column: **TRB-624**, P/N TR-603075


Dimensions: 75m x 0.53mm x 3.0  $\mu$ m

Injection: mixture of solvents (wet needle), split 1:100, 250°C

Carrier gas: He, constant pressure 8 psi (55.7 Kpa)

Oven temperature: 40°C(15 min) @ 15°C/min to 75°C(15 min)

Detector: FID, 250°C

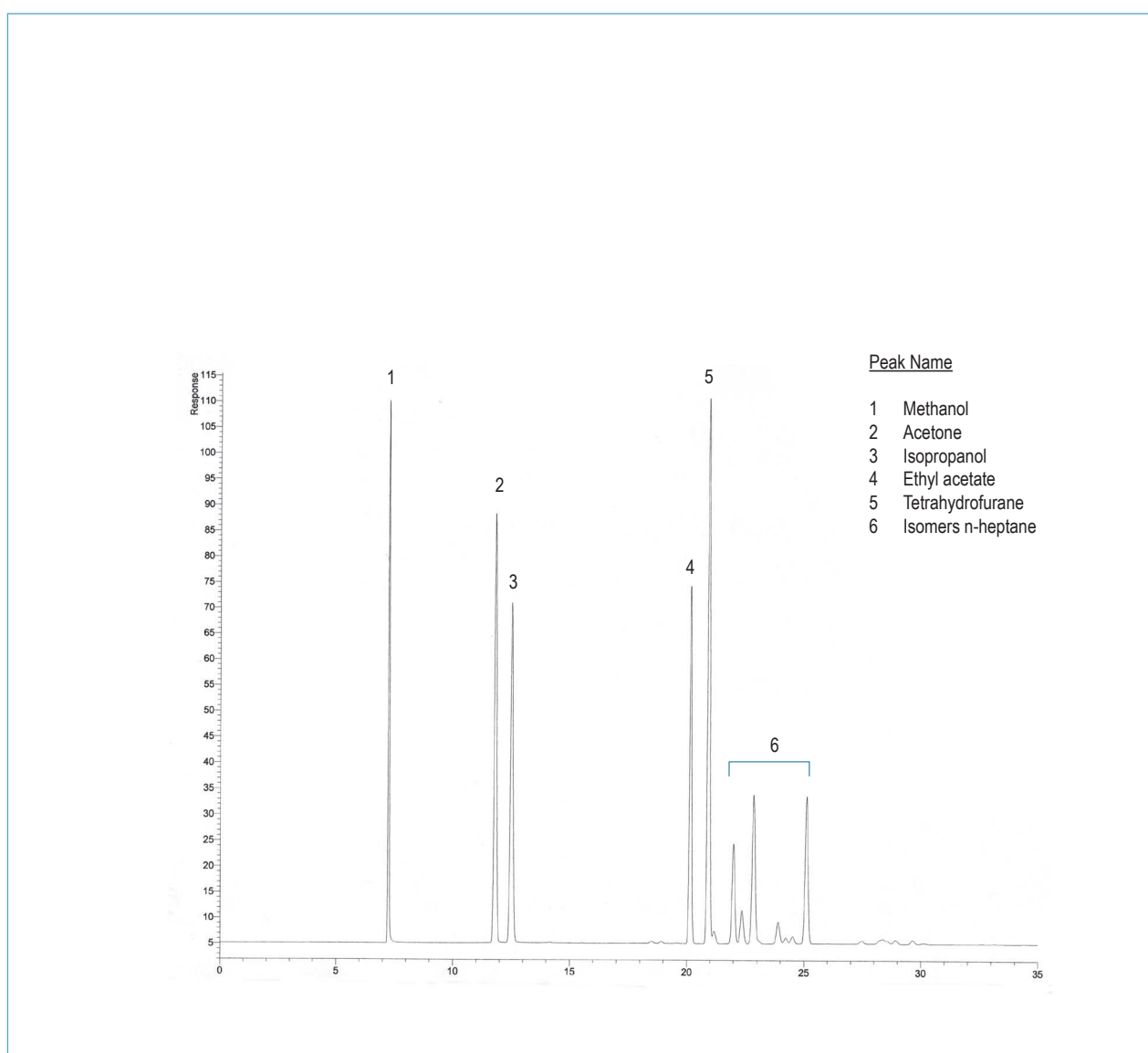

TKG 1201

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: [comercial@teknokroma.es](mailto:comercial@teknokroma.es)