

METHANOL AND HIGHER ALCOHOLS

Column: **TRB-624**, P/N TR-601432

Size: 30m x 0.25mm x 1.4 μ m

Injection: 1 μ L split 1:50, 250°C



Carrier gas: He, 1mL/min

Program temperature: 40°C (5min) @ 20°C/min to 200°C (3min)

Detector: FID, 260°C

Chromatogram provided by Joan Garcia from INCAVI (Vilafranca del Penedés, Barcelona)

RT (min)	Peak Name
3.21	Methanol
4.66	Ethanol
7.29	1-Propanol
8.16	2-Butanol
8.86	Isobutanol
9.60	1-Butanol
10.87	Isoamylic alcohol
10.92	Isoamylic alcohol

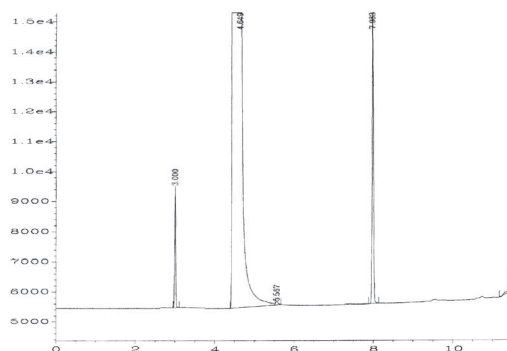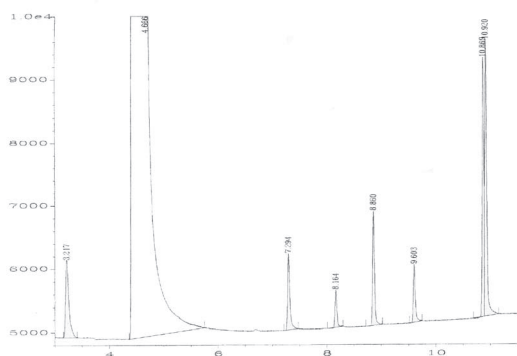

TKG 1242

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es