

REACTION PRODUCTS OF 1,1-DIETHOXY BUTANE

Column: **TRB-624**, P/N TR-603035

Size: 30m x 0.53mm x 3.0µm


Injection: 0.5 µl, split 1:5, 260°C

Carrier Gas: He, 6 psi

Program temperature: 40°C (6min) @ 30°C/min to 200°C (5min)

Detector: FID, 260°C

Chromatogram provided by Ion Aguirre from Escuela Superior de Ingeniería de Bilbao (Spain)

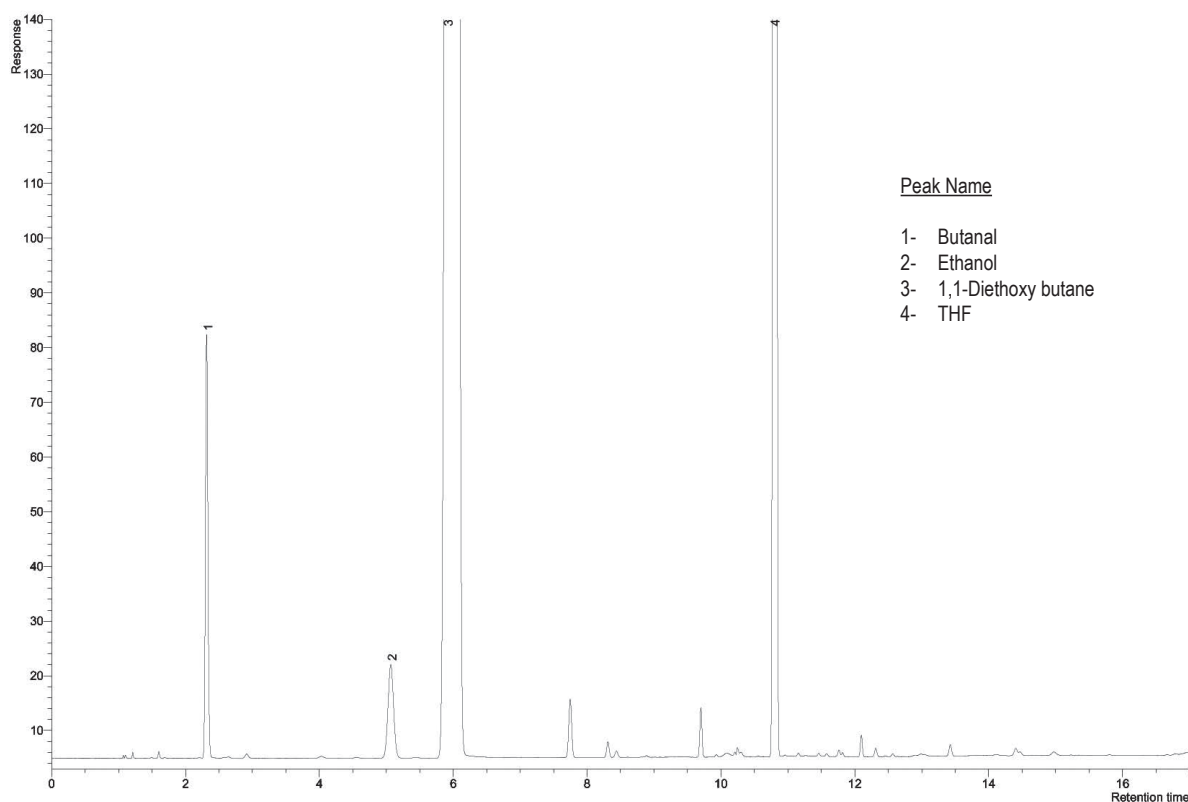

TKG 1238

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es