

TRB-14

Column: TRB-14, P/N TR-801263


Dimensions: 60m x 0.32mm x 1.2µm

Injection: 0.1µL Amines mixture (0.02-0.2% in Acetonitrile), split 1:25, 200°C

Carrier gas: He, constant pressure 14 psi (96.5 Kpa)

Oven temperature: 40°C(2min) @ 10°C/min to 250°C

Detector: FID, 280°C

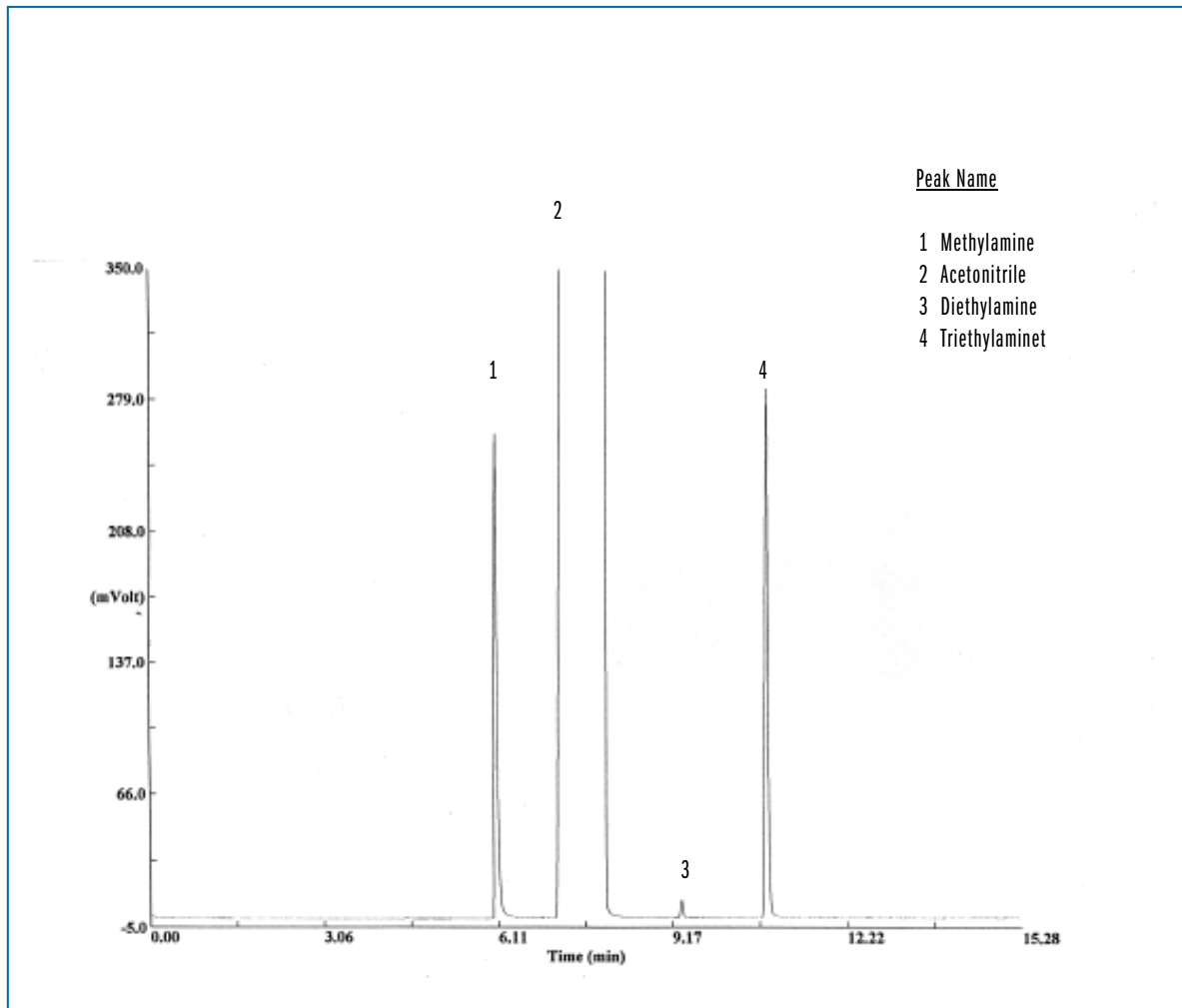

TKG 1123

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain

Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es