

ANALYSIS OF SIMAZIN

Column: **TRB-50**, P/N TR-500232


Dimensions: 30m x 0.25mm x 0.25 μ m

Injection: 0.2 μ L standard (20 ng on column), split 1:10, 280°C

Carrier gas: H₂, constant pressure 11 psi (75.8 KPa)

Oven temperature: 200°C (Isothermal)

Detector: FID, 280°C

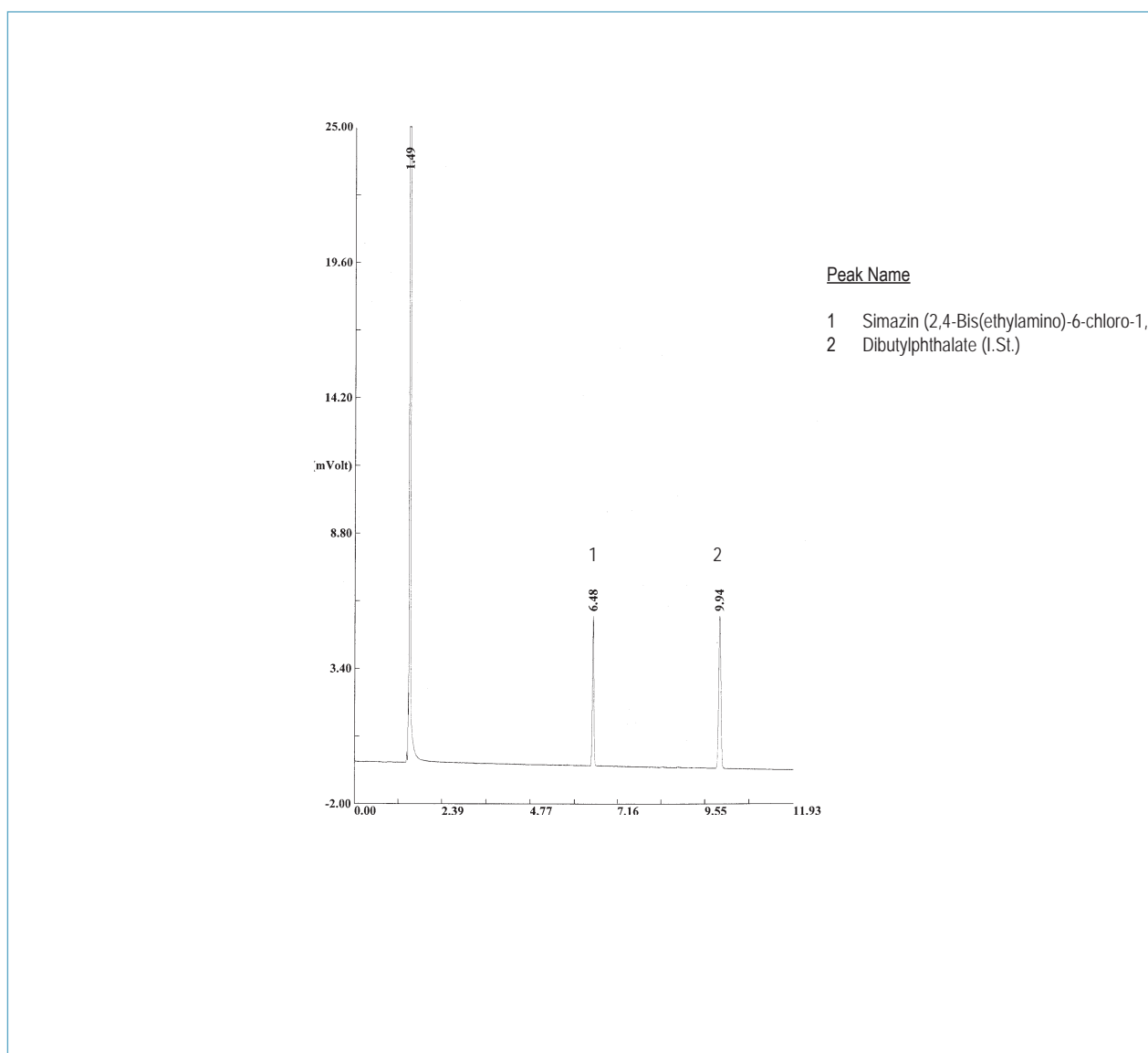

TKG 1001

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es