


**TRB-50ht**

Column: **TRB-50ht**, P/N TR-531312  
 Dimensions: 15m x 0.25mm x 0.15mm  
 Injection: 0.2mL Triglycerides Palm Oil in Isooctane(50 mg/mL), split 1:12  
 Carrier gas: H2, constant pressure, 9psi (56 KPa)  
 Oven temperature: 340°C(1min)@0.5°C/min to 355°C(5min)  
 Injector: 380°C (high temperature septum)  
 Detector: FID, 380°C

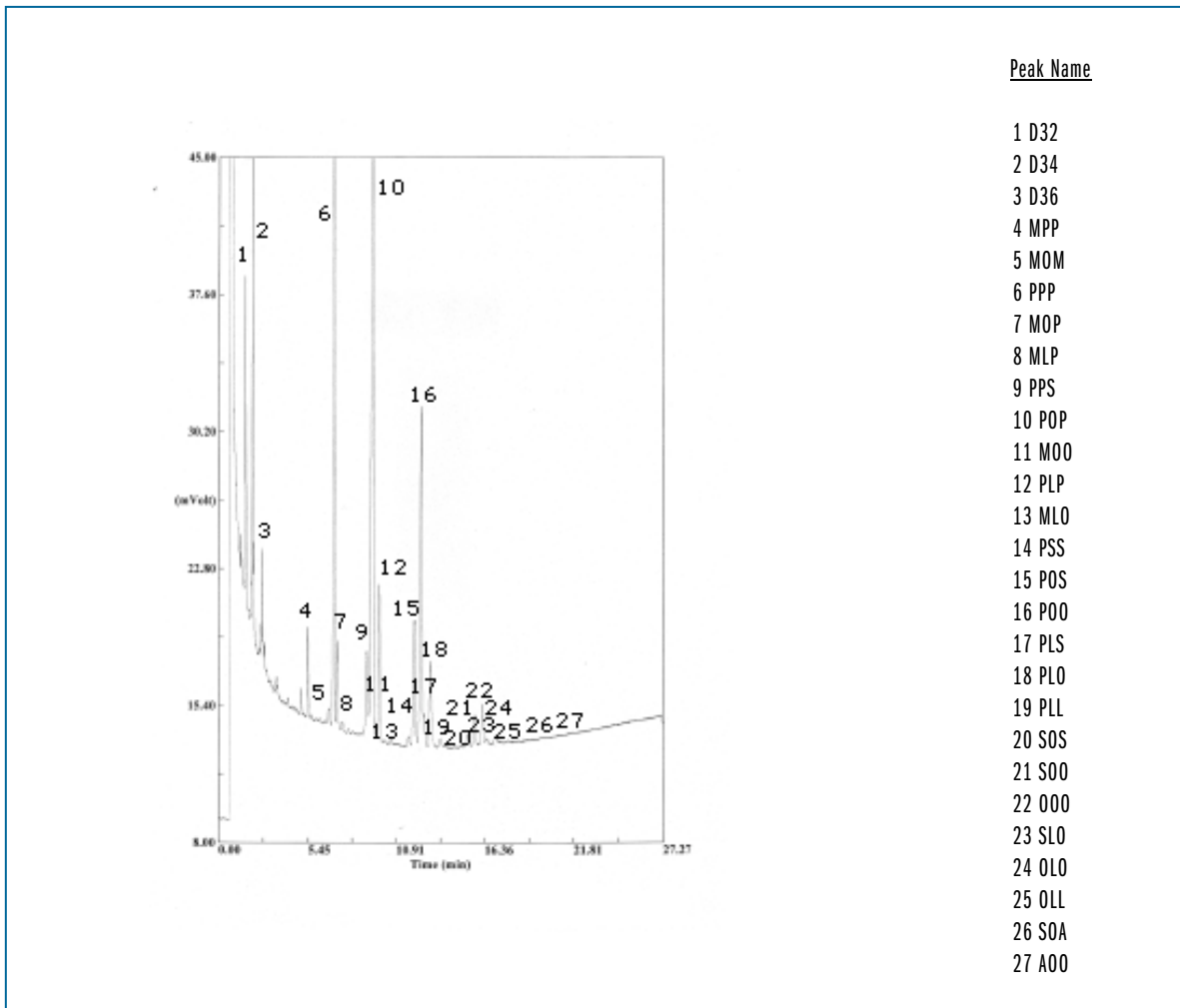

TKG 1128