

AMINOACIDS

Column: **TRB-50ht**, P/N TR-531332

Size: 30m x 0.25mm x 0.15µm


Injection: 2 µl standard AA-S-18 Sigma (2.5µmol/ml), split, 300°C

Carrier gas: He, 1mL/min

Program temperature: 50°C (2min) @ 30°C/min to 350°C (3min)

Detector: MS Polaris Q, EI, 200°C, transfer line 200°C

Chromatogram provided by Antonio Tintó from Moehs S.A.

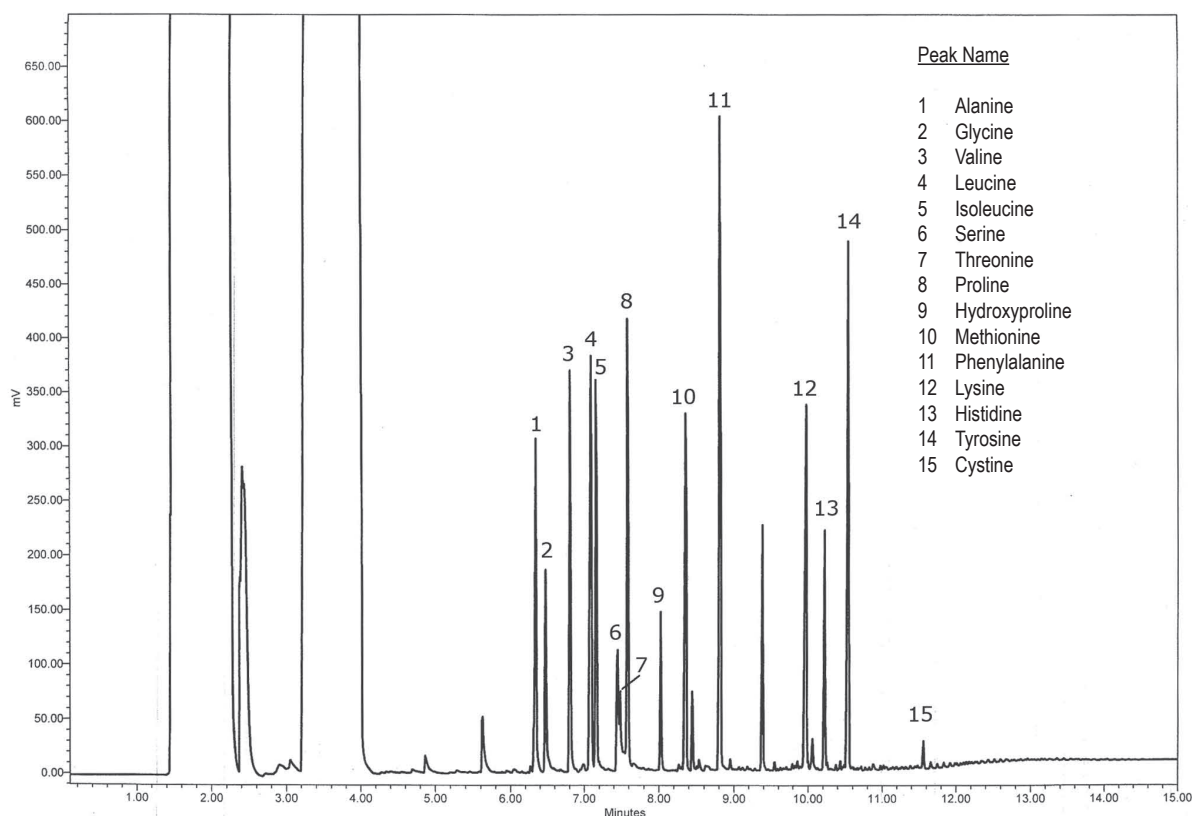

TKG 1215

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es