


TRB-PAG

Column: TRB-PAG, P/N TR-550232
Dimensions: 30m x 0.25mm x 0.25µm
Injection: 1 µL Test Grob, split 1:25, 260°C
Carrier gas: H2, constant pressure 11 psi(75.8 Kpa)
Oven temperature: 40°C @ 6°C/min to 230°C(5min)
Detector: FID, 260°C

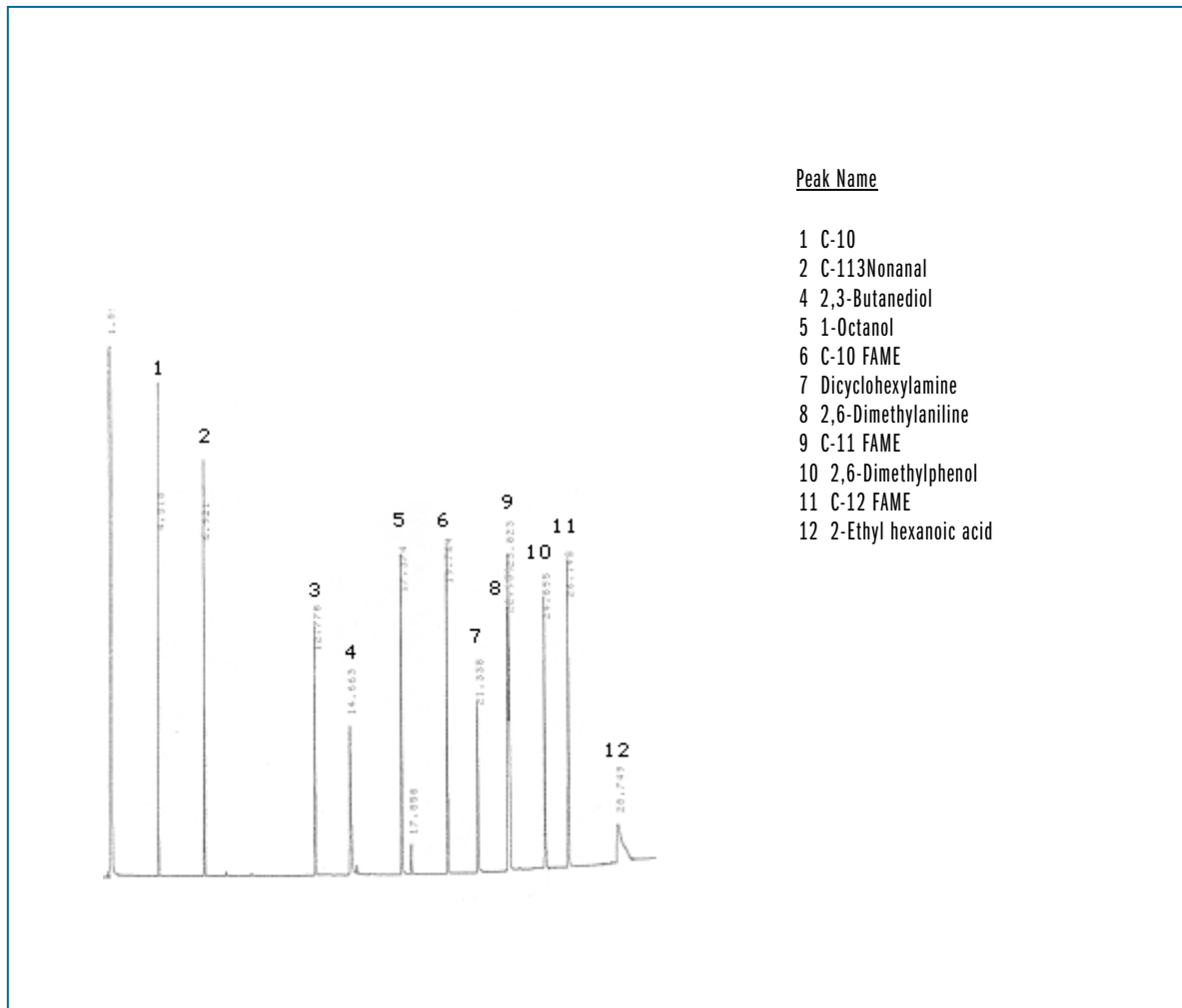

TKG 1130

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es