


SEPARATION OF SOLVENTS

Column: **TRB-WAX**, P/N TR-141253
 Dimensions: 50m x 0.32mm x 1.2µm
 Injection: 1 µL standard (500 ng/mL comp.), split 1:25, 260°C
 Carrier gas: He, constant pressure 12 psi (82.7 Kpa)
 Oven temperature: 65°C(7min) @ 4°C/min to 117°C
 Detector: FID, 260°C

*Chromatogram provided by
Jaume Piedrabuena from Danisco*

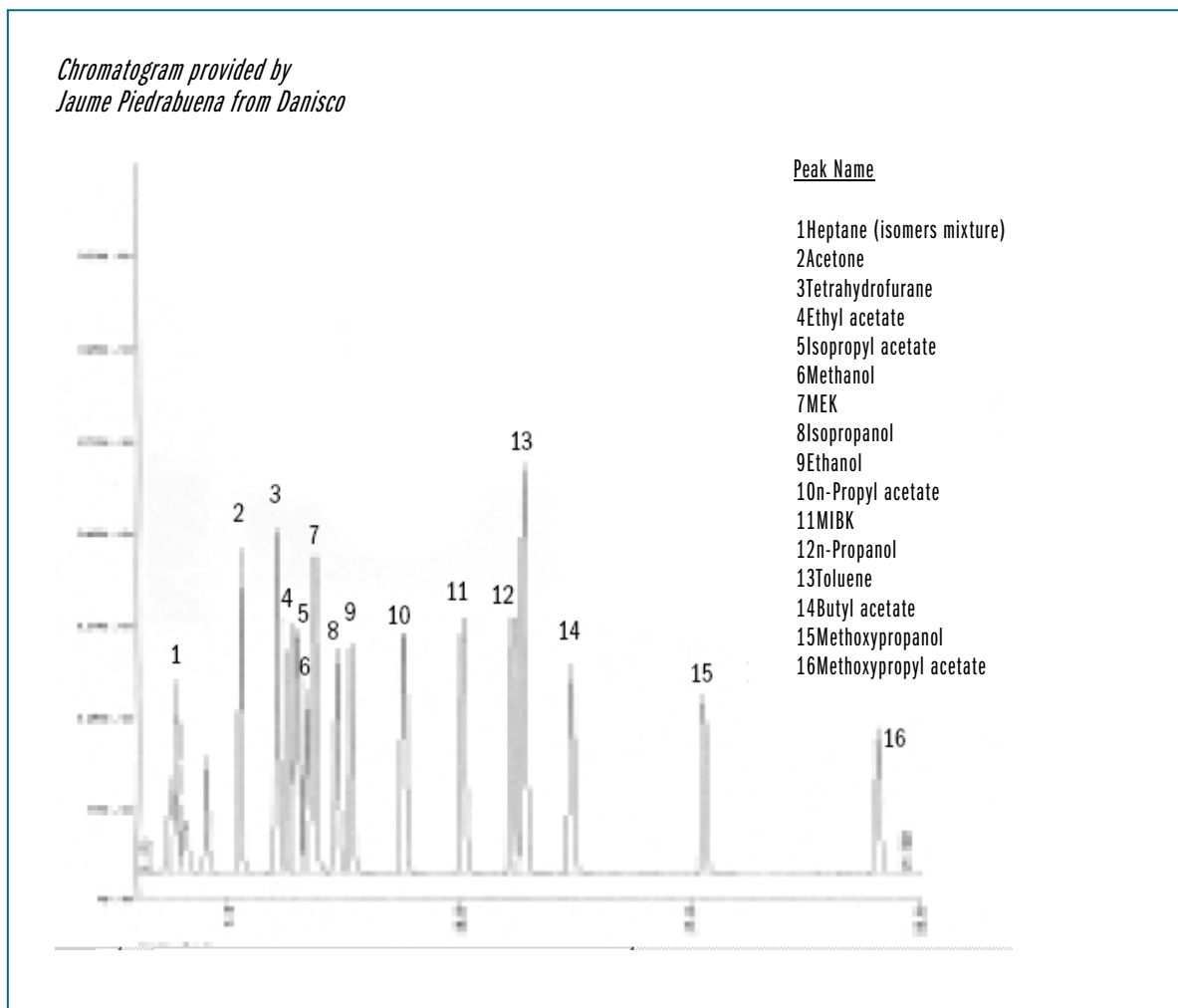

TKG 1052