


## ANALYSIS OF VOLATILES IN OLIVE OIL

Column: **TRB-WAX**, P/N TR-140262  
Dimensions: 60m x 0.25mm x 0.25µm  
Injection: 1µL Olive oil, SPME (SP-57329-U), 260°C  
Carrier gas: H<sub>2</sub>, 1.5mL/min  
Oven temperature: 40°C @ 2°C/min to 200°C  
Detector: FID, 280°C

*Chromatogram provided by Departamento de Caracterización y Calidad de los Alimentos, CSIC (Sevilla)*

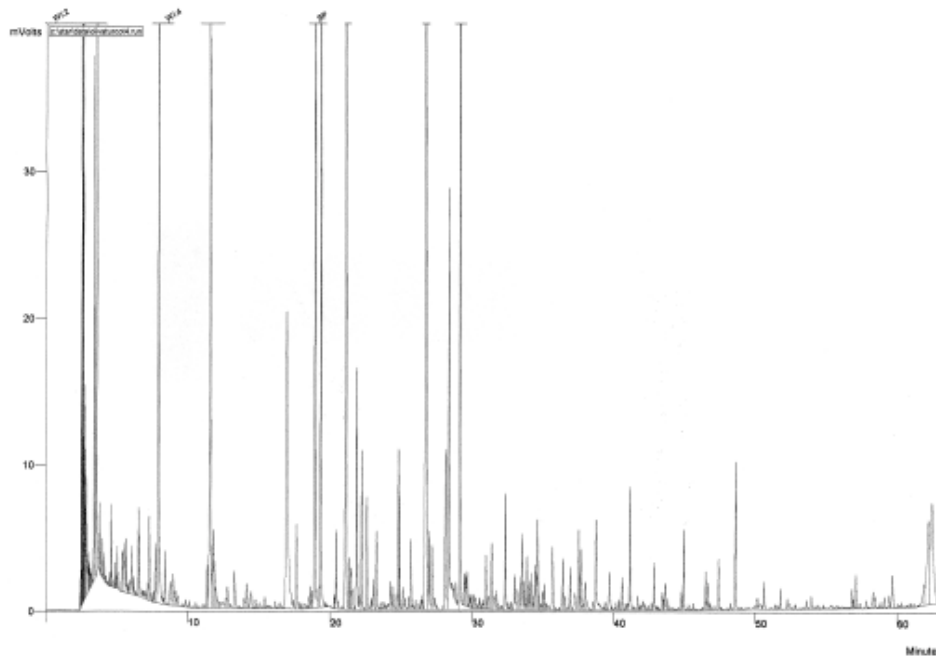

TKG 1091

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es