

**TRB-WAX**

Column: TRB-WAX, P/N TR-140141


Dimensions: 10m x 0.10mm x 0.10µm

Injection: 1µl standard FAMES, (200ng/comp), split 1:50, 280°C

Carrier gas: H<sub>2</sub>, constant pressure 50 psi (344.5 Kpa)

Oven temperature: 40°C(1min) @ 25°C/min to 195°C @ 3°C/min to 205°C @ 8°C/min to 230°C(1min)

Detector: FID, 280°C

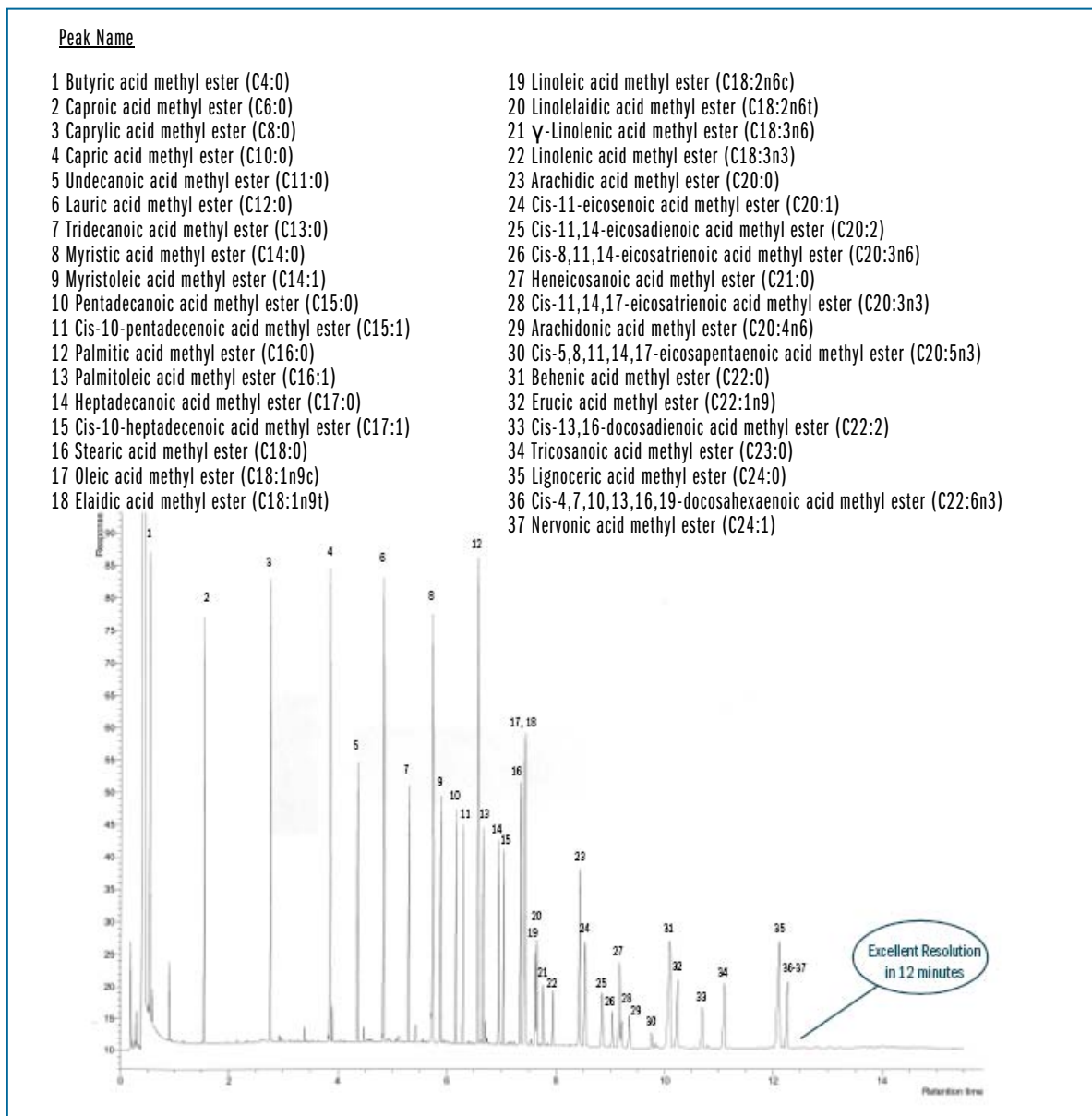

TKG 1145