


## LAVANDER FLAVOR

Column: **TRB-WAX**, P/N TR-140232  
Dimensions: 30m x 0.25mm x 0.25µm  
Injection: 1 µL, split  
Carrier gas: He, 90 KPa  
Oven temperature: 80°C @ 4°C/min to 230°C(20 min)  
Detector: FID, 260°C

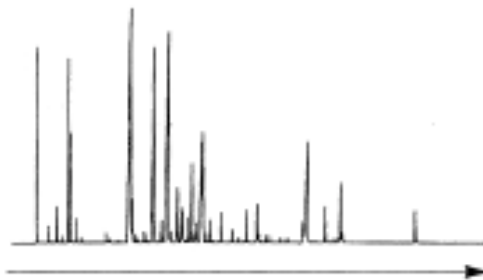

*Chromatogram provided by C. Ibañez from Lueta, S.A, Barcelona.*

TKG 1181

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es