


FLAVORS (LAVANDER, ESSENTIAL OIL)

Column: **TRB-WAX**, P/N TR-140232
Dimensions: 30m x 0.25mm x 0.25 μ m
Injection: 1 μ L, split
Carrier gas: He, 90 KPa
Oven temperature: 80°C @ 4°C/min to 230°C(20 min)
Detector: FID, 260°C

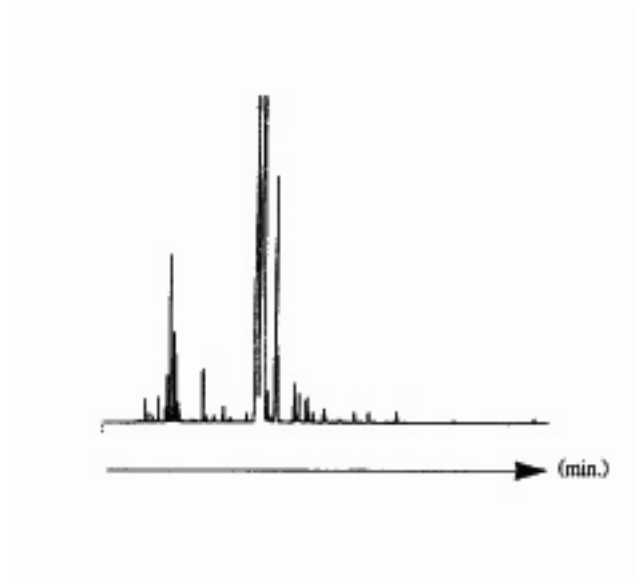

Chromatogram provided by C. Ibañez from Lucta, S.A, Barcelona.

TKG 1182

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es