

Column: **TRB-WAX**, P/N TR-140232

Size: 30m x 0.25mm x 0.25µm

Injection: split 1:30, 250°C (liner SPME, fibra: 2cm 50/30 µm DVB/Carboxen/PDMS)

Sample: 5 µl soap in 4 mL of water (26.6% NaCl)


Carrier gas: He, 1mL/min

Program temperature: 32°C (5min) @ 4°C/min to 250°C (5min)

Detector: MS KONIK-TECH, Mode EI+ (70 eV), range mass 35-300, Scan time 35ms,

Source Temperature 120°C, Interface Temperature 250°C, Photomultiplier 750V.

Chromatogram provided by Nieves Sarrión de KONIK-TECH, S.A

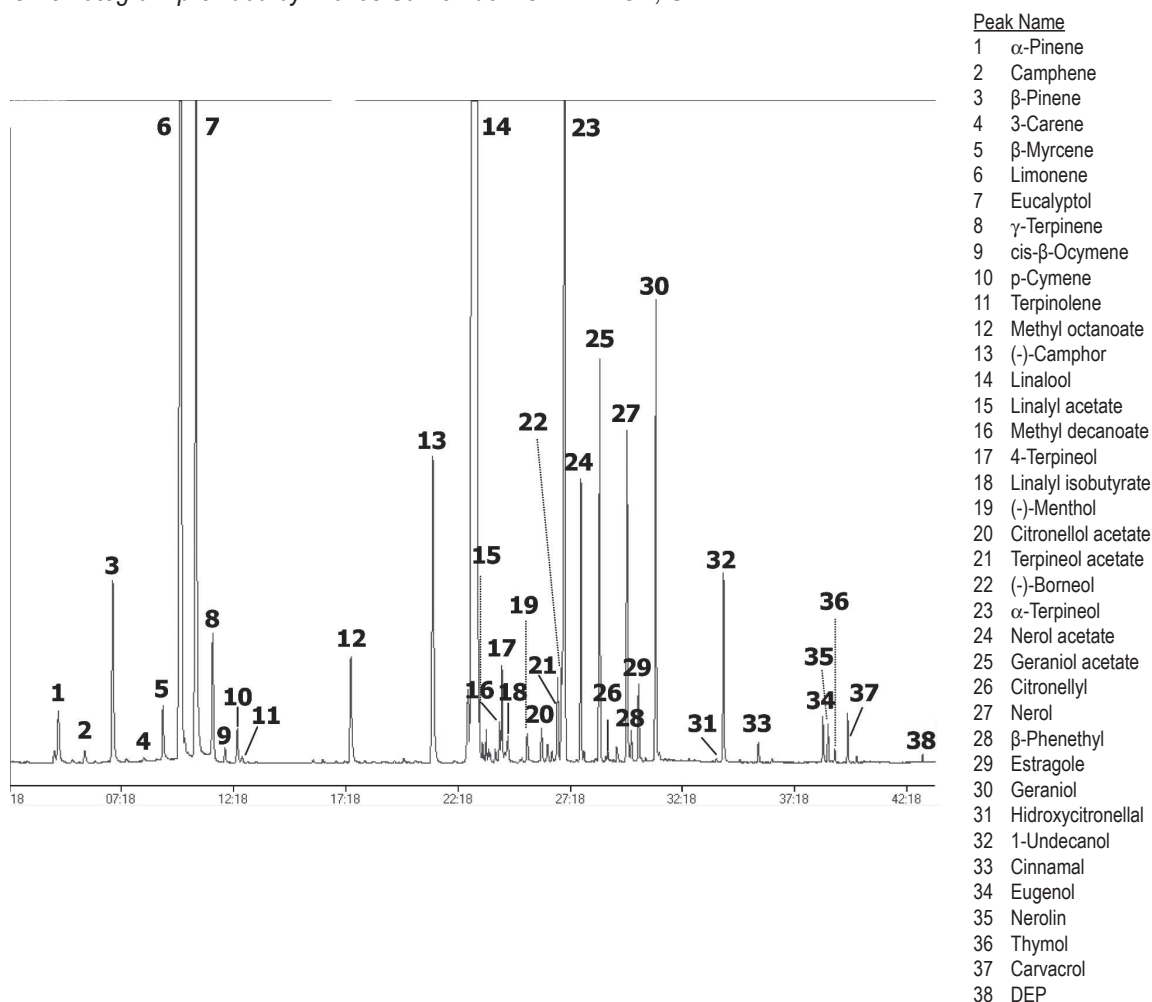

TKG 1228

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es