

SEPARATION OF ALDEHYDES

Column: **Meta .WAX**, P/N TR-810532


Dimensions: 30m x 0.25mm x 0.5 μ m

Injection: wet needle (neat), split 1:100, 250°C

Carrier gas: He, 12 psi (82.7 KPa)

Oven temperature: 40°C @ 10°C/min to 120°C

Detector: FID, 260°C

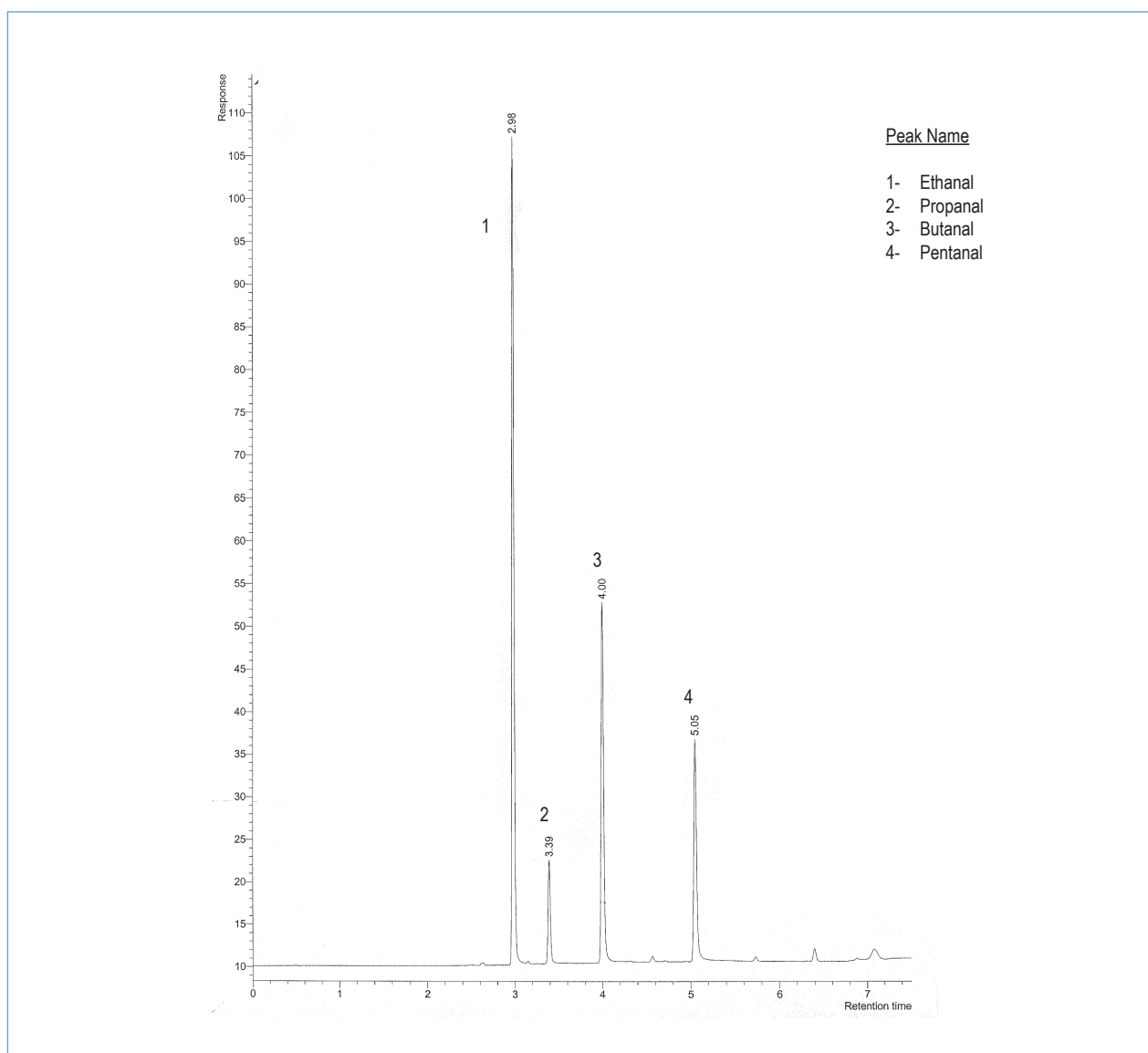

TKG 1018

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es