


ANALYSIS BERGAMOL

Column: Meta.WAX, P/N TR-810532
Dimensions: 30m x 0.25mm x 0.5µm
Injection: 1 µL standard 0.3% v/v Bergamol/Diethyl Phtalate in Ethanol, split 1:50, 260°C
Carrier gas: H₂, 12 psi (82.7 KPa).
Oven temperature: 35°C(10min) @ 8°C/min to 220°C(20min)
Detector: FID, 260°C

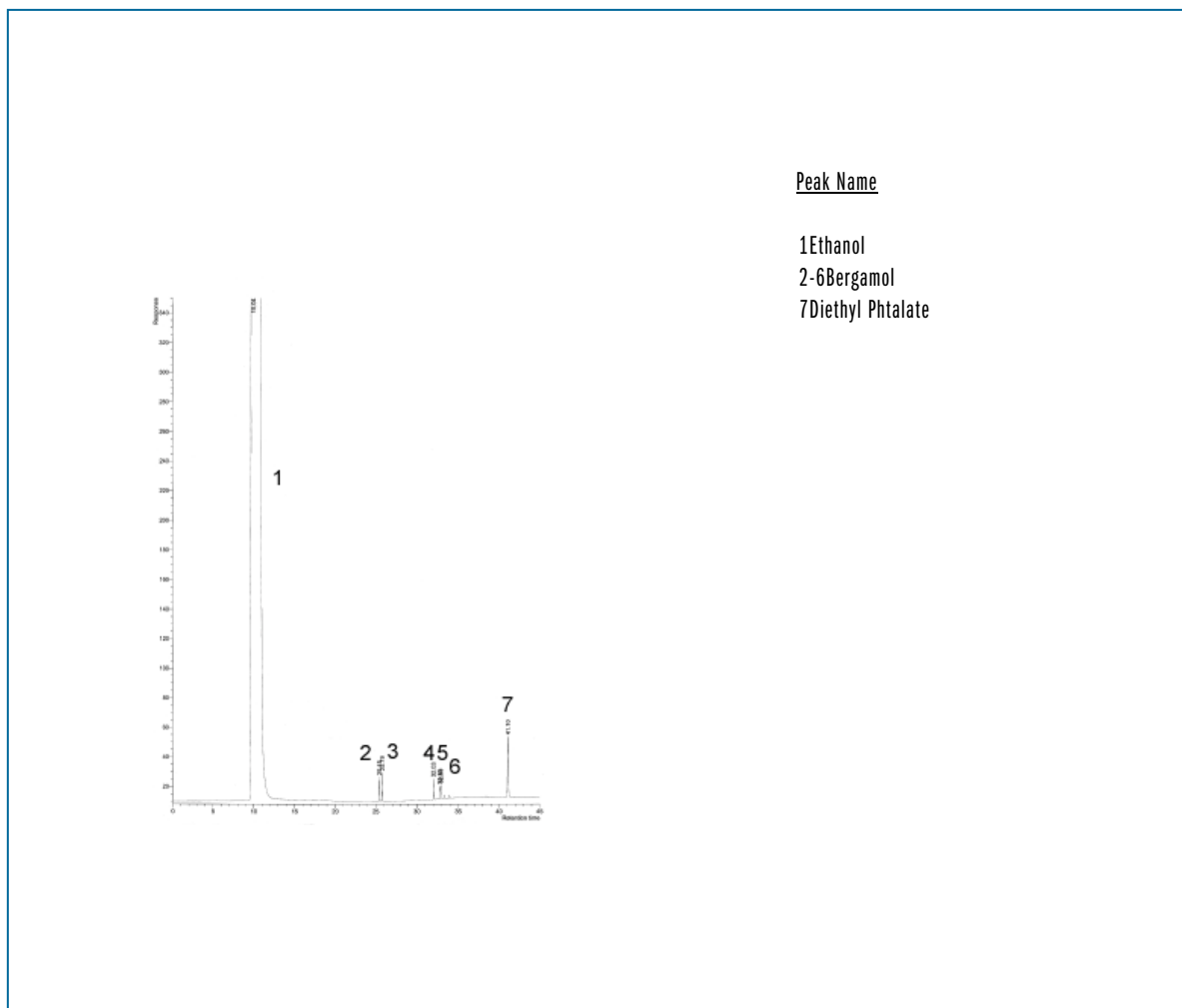

TKG 1050

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es