


Meta.WAX

Column: Meta.WAX, P/N TR-810532
 Dimensions: 30m x 0.25mm x 0.5µm
 Injection: 1µL Standard Alcohols, split 1:25, 260°C
 Carrier gas: He, 12 psi (82.7 KPa)
 Oven temperature: 35°C(10min) @ 8°C/min to 220°C (15min)
 Detector: FID, 260°C

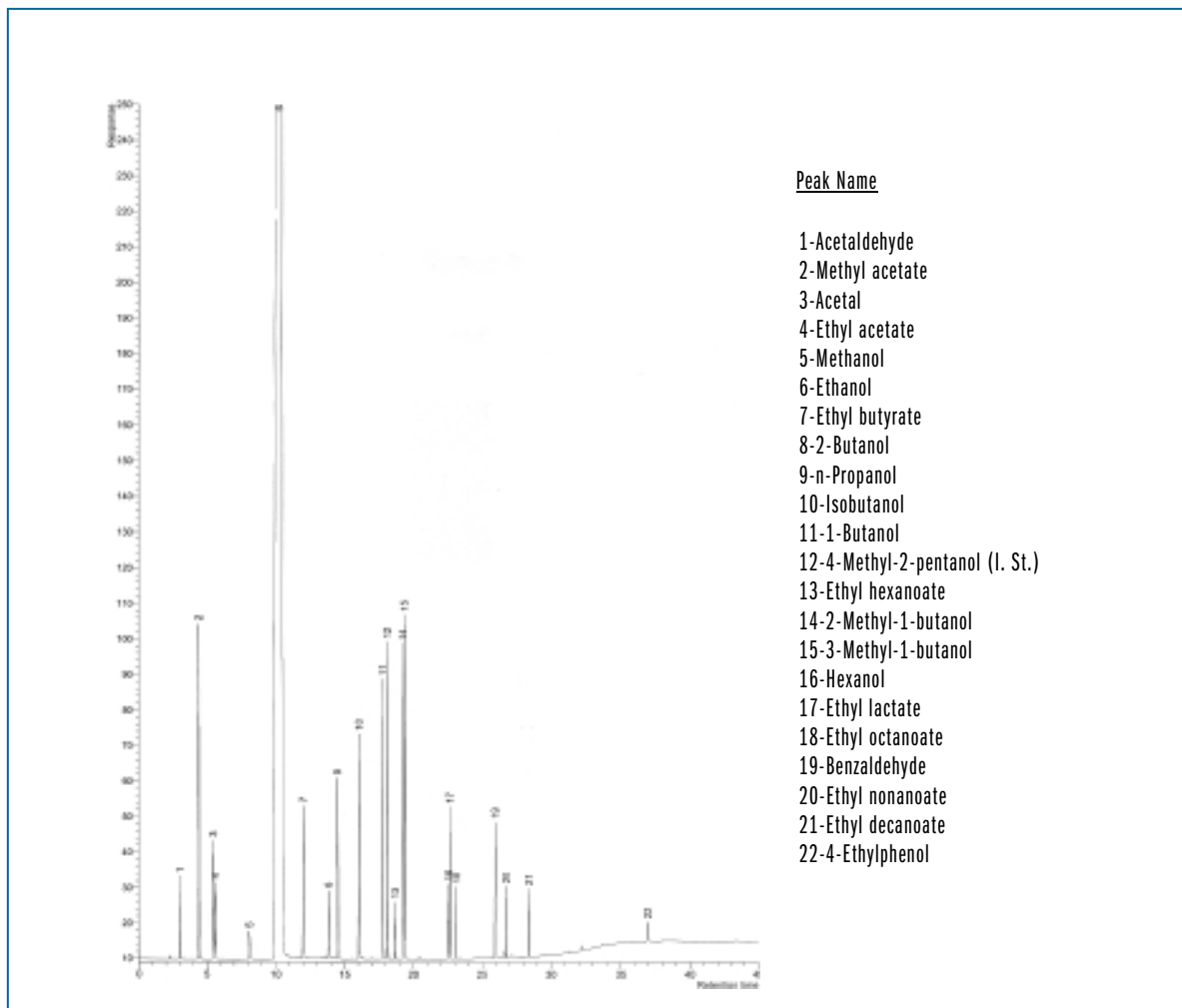

TKG 1135

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es