


SEPARATION OF SOLVENTS

Column: **TR-Meta.WAX 400**, P/N TR-402153
Dimensions: 50m x 0.32mm x 0.20µm
Injection: 0.5µL standard(split 1:50), 200°C
Carrier gas: He, 16 psi (110.2 Kpa)
Oven Temperature: 30°C(2') @1°C/min to 80°C
Detector Temperature: FID, 200°C

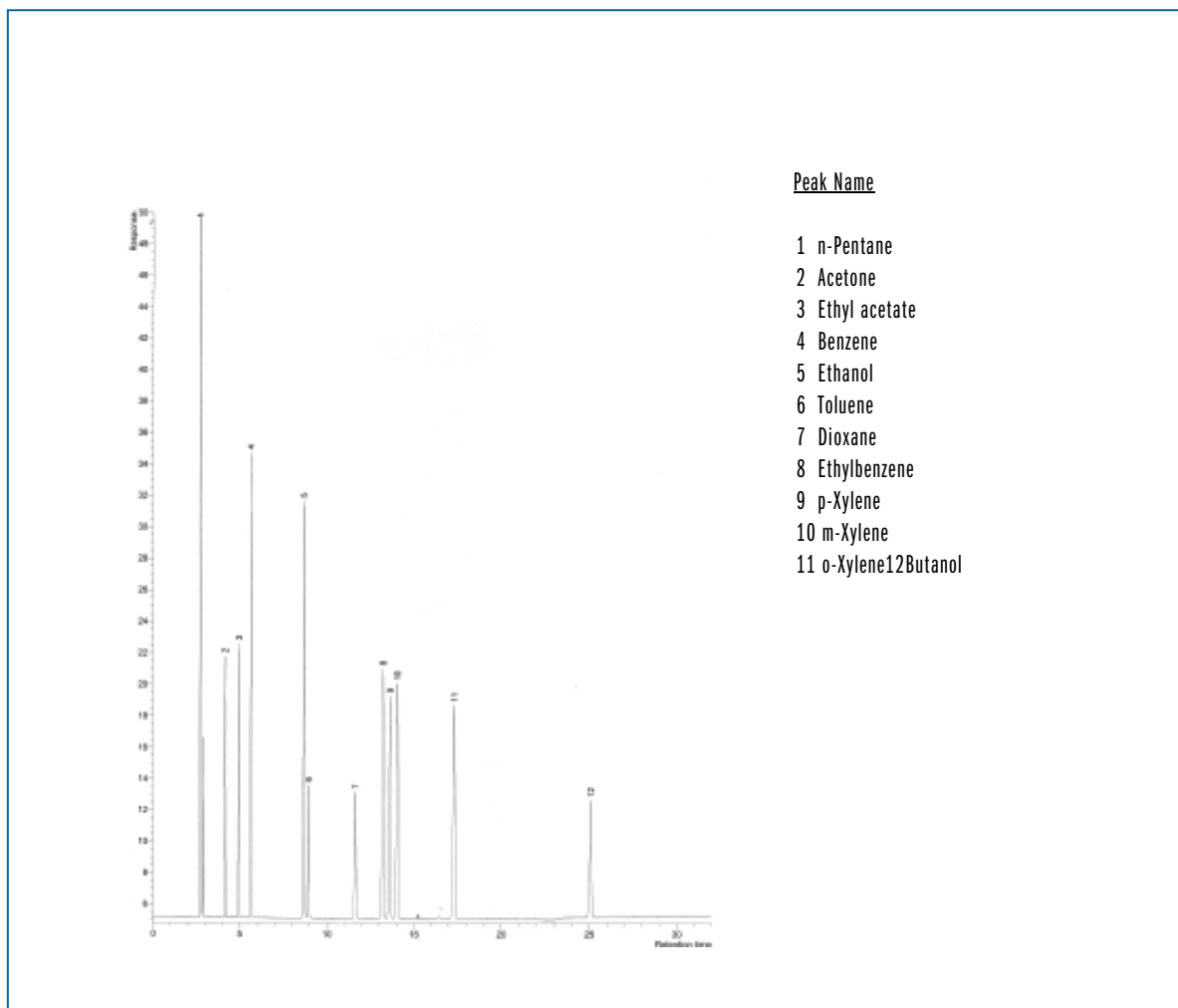

TKG 1060

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es