

## SEPARATION OF METHYL ESTERS IN SOYA OIL

Column: **TR-CN100**, P/N TR-882162


Size: 60m x 0.25mm x 0.20µm

Injection: 1µL Total FAMES en CH<sub>2</sub>Cl<sub>2</sub> (30 mg/mL), split 1:100, 280°C

Carrier gas: He 25 psi

Program temperature: 185°C

Detector: FID, 280°C

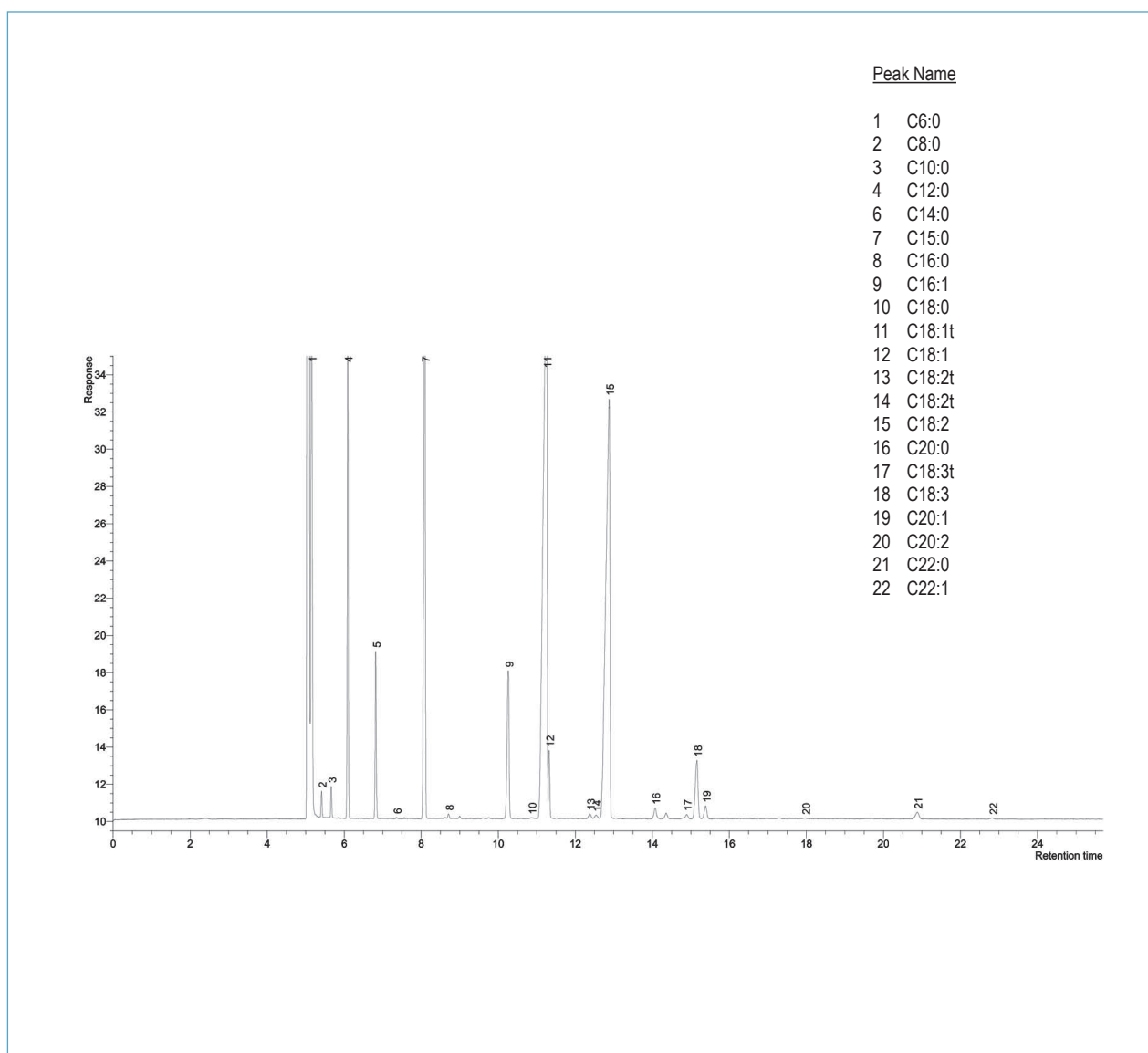

TKG 1239

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es