

37 FAME MIX

Column: **TR-CN100**, P/N TR-8816H4


Size: 75m x 0.18mm x 0.14µm

Injection: 1µL 37 FAME MIX (47885-U), split 1:100, 250°C

Carrier gas: H2 40 cm/s (140°C)

Program temperature: 140°C(5min) @ 4°C/min to 240°C(2min)

Detector: FID, 250°C

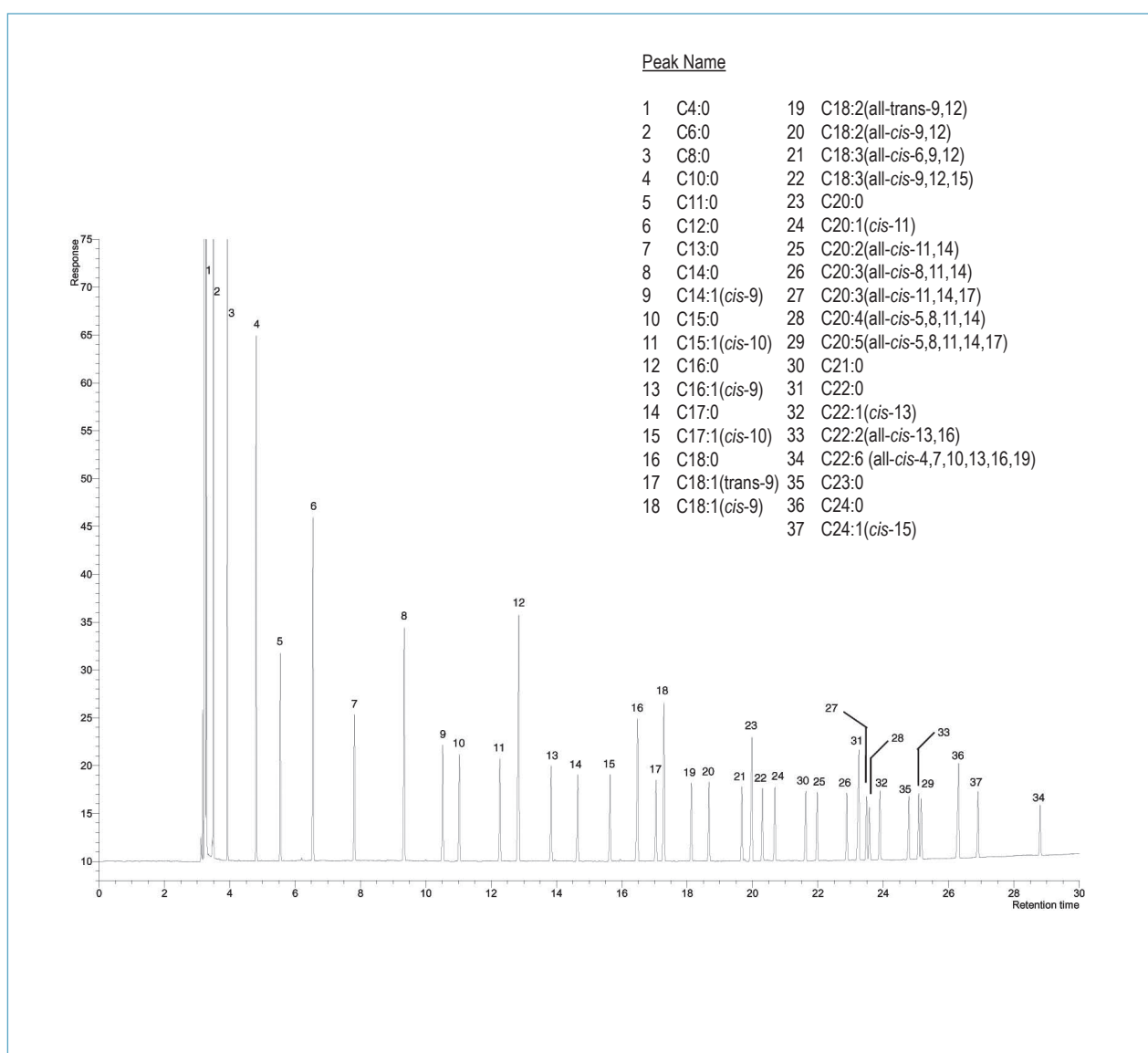

TKG 1250

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es