


TR-CRESOL

Column: **TR-CRESOL**, P/N TR-712162
 Dimensions: 60m x 0.25mm x 0.20µm
 Injection: 1 µL standard Cresols (5000 ng/mL comp.), split 1:25, 150°C
 Carrier gas: H₂, constant pressure 24 psi (165 Kpa)
 Oven temperature: 130°C
 Detector: FID, 150°C

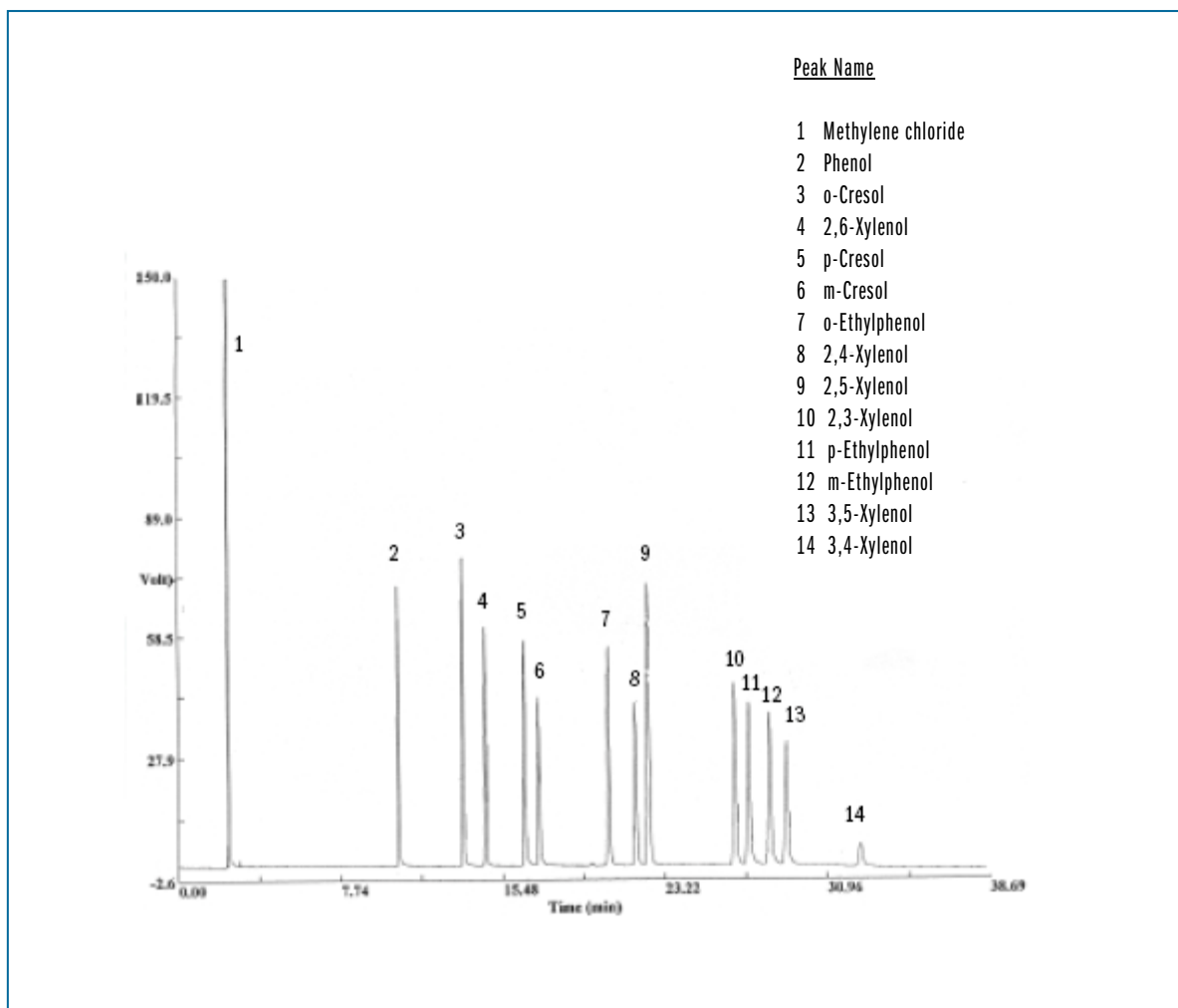

TKG 1137

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es