


TR-17

Column: TR-17 P/N TR-502045
Dimensions: 10m X 0,53 mm X 2.0 µm
Injector: 260°C
Carrier gas: He, 6.5 psi
Injection: 1ml standard, split (1:4)
Oven Temperature: 180°C
Detector: FID, 280°C

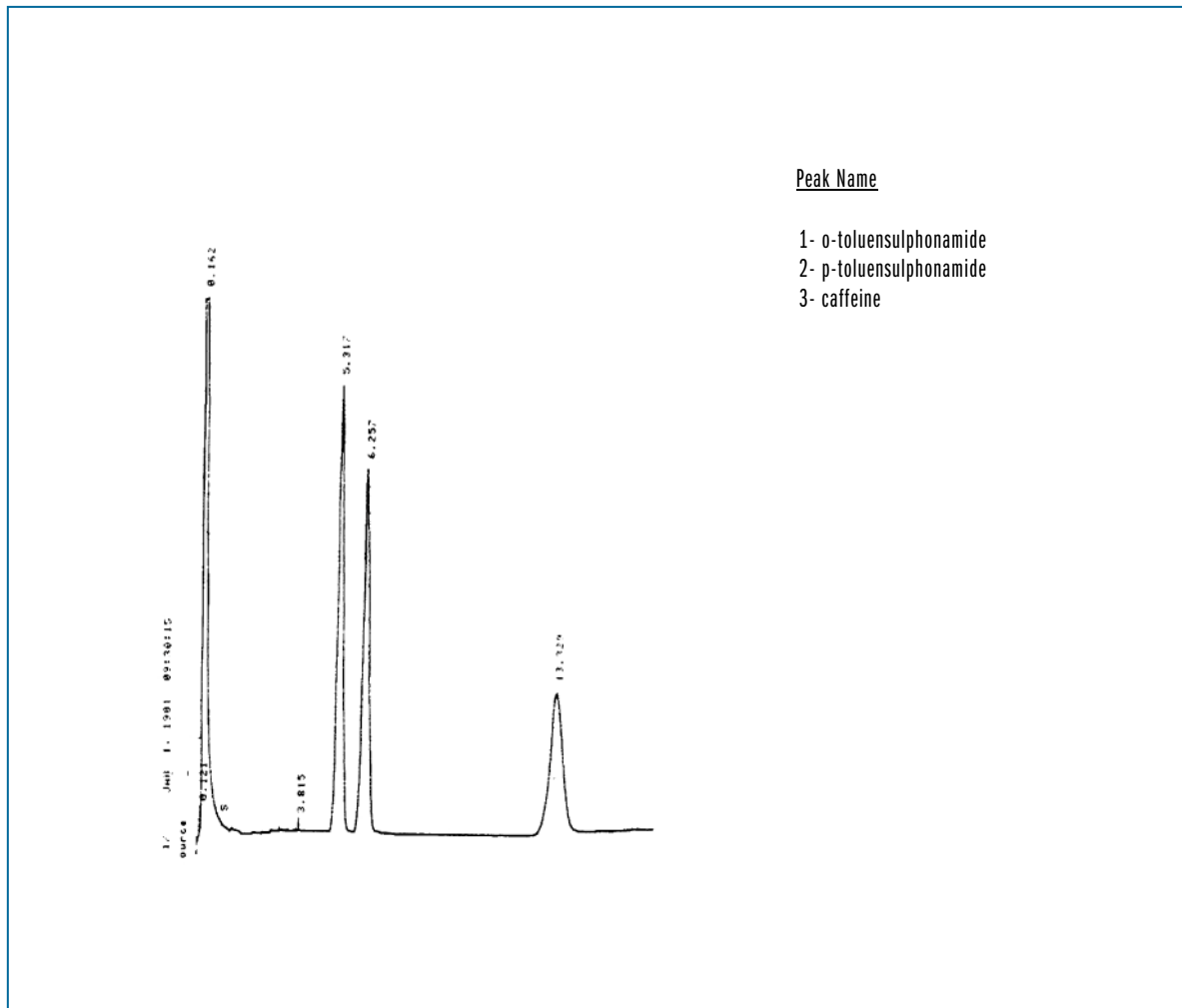

TKG 1138

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es