


TRB 608

Column: **TRB-608**, P/N TR-360232
 Dimensions: 30m x 0.25mm x 0.25mm
 Injection: 1 µL Grob Test, split 1:100, 280°C
 Carrier gas: H₂, constant pressure 11 psi (75.8 Kpa)
 Oven temperature: 40°C @ 6°C/min to 300°C(5min)
 Detector: FID, 340°C

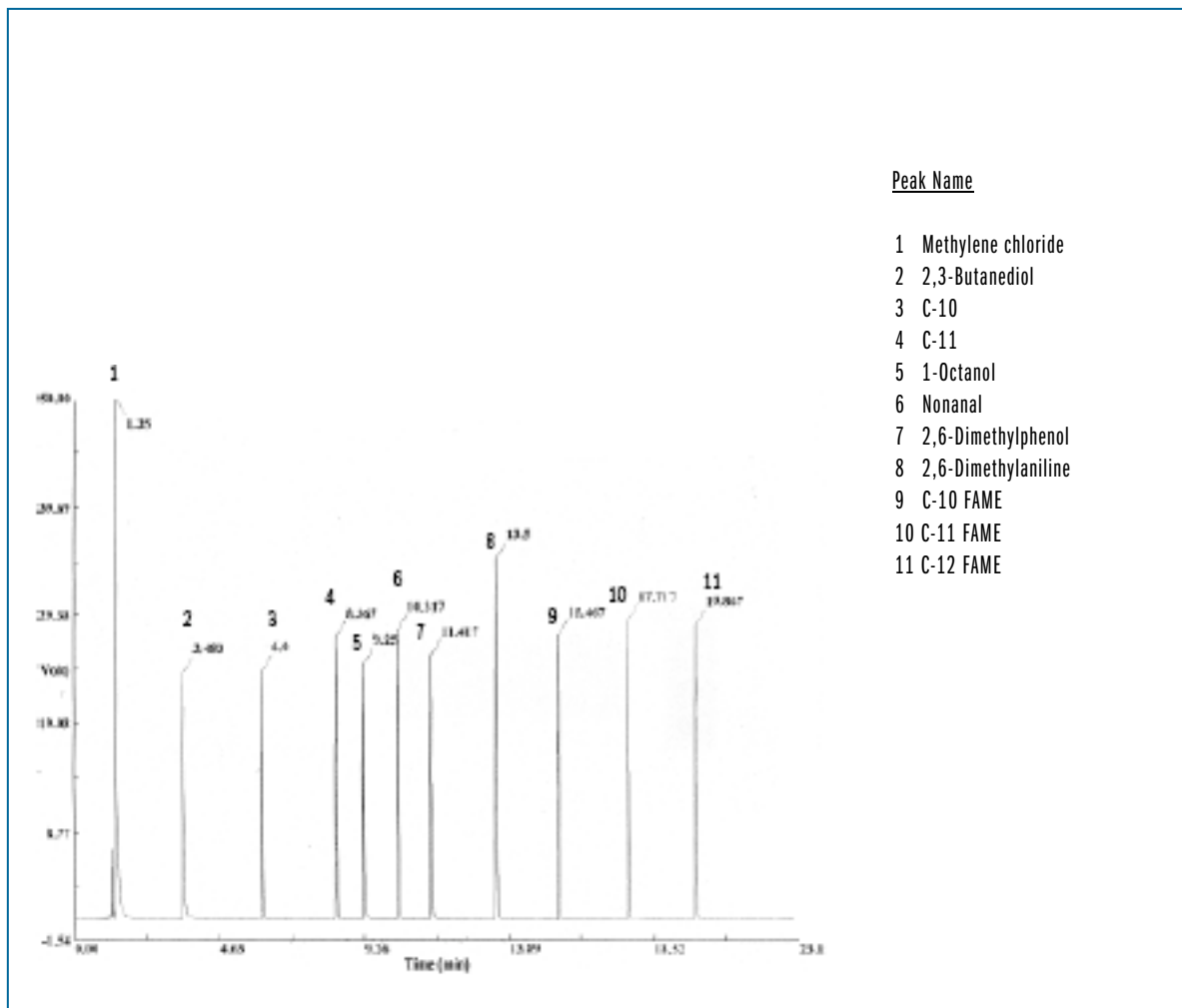

Peak Name

- 1 Methylene chloride
- 2 2,3-Butanediol
- 3 C-10
- 4 C-11
- 5 1-Octanol
- 6 Nonanal
- 7 2,6-Dimethylphenol
- 8 2,6-Dimethylaniline
- 9 C-10 FAME
- 10 C-11 FAME
- 11 C-12 FAME

TKG 1140

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es