

ANALYSIS OF PESTICIDES

Column: **TRB-608**, P/N TR-360232


Size: 30m x 0.25mm x 0.25 μ m

Carrier Gas: He, 18 psi (124 KPa), 35 cm/s

Program Temperature: 150°C(4 min) @ 8°C/min to 290°C(10 min)

Injection: 1 μ L Test Pesticides (40 a 80 μ g/l), splitless, 175°C

Detector: ECD, 300°C

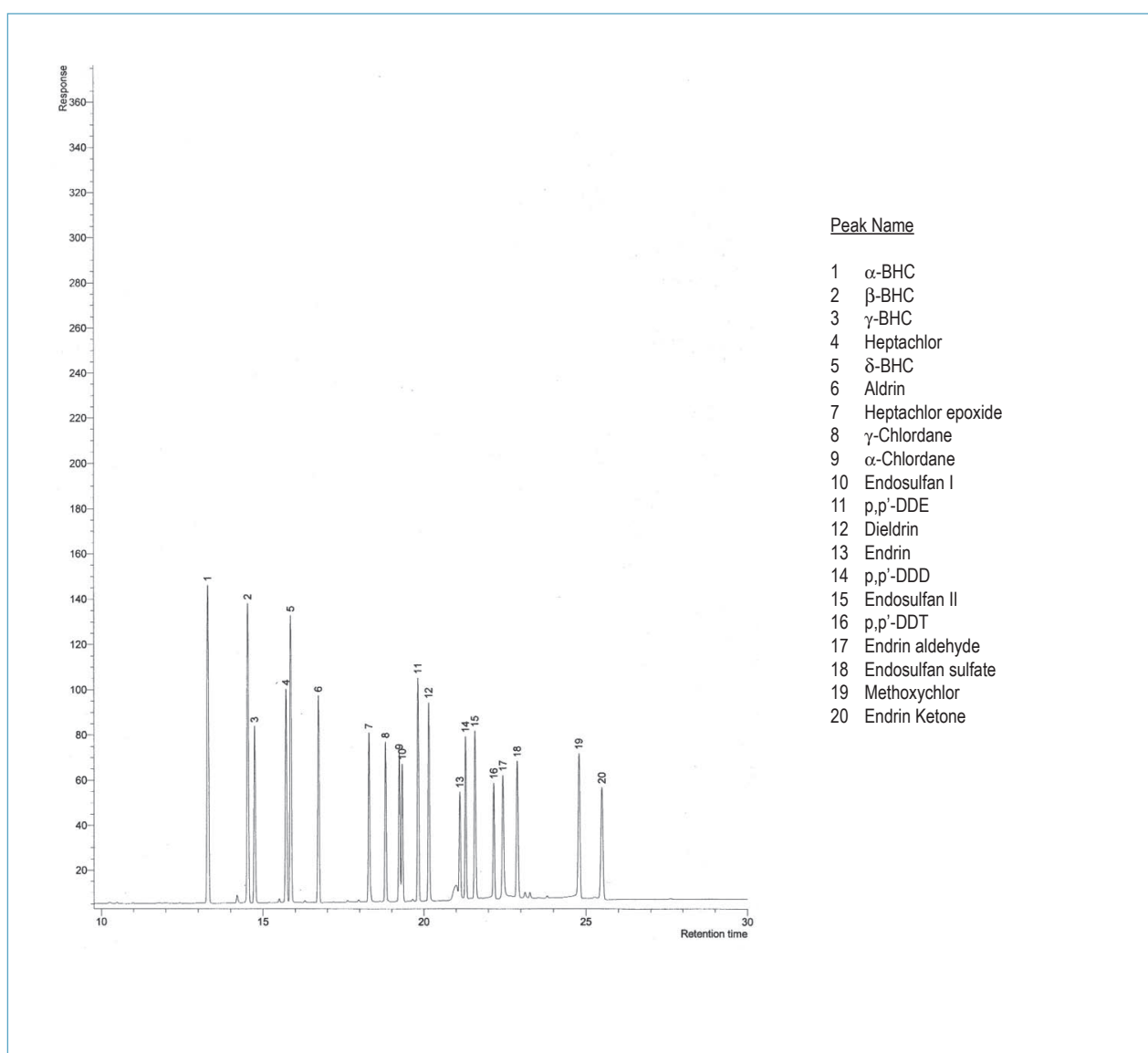

TKG 1211

www.teknokroma.com

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es