

SAPIENS-624MS: Residual Solvents

Column: SAPIENS-624MS P/N TR-591833


Dimensions: 30 m x 0,32 mm x 1,8 µm

Injection: 0.1 µl Residual Solvent mix (0,1%) 260°C, Split: 1:50

Carrier gas: He, 10 psi

Oven temperature: 35°C (20 min) to 150°C (2min) @ 12°C/min

Detector: FID @ 260°C

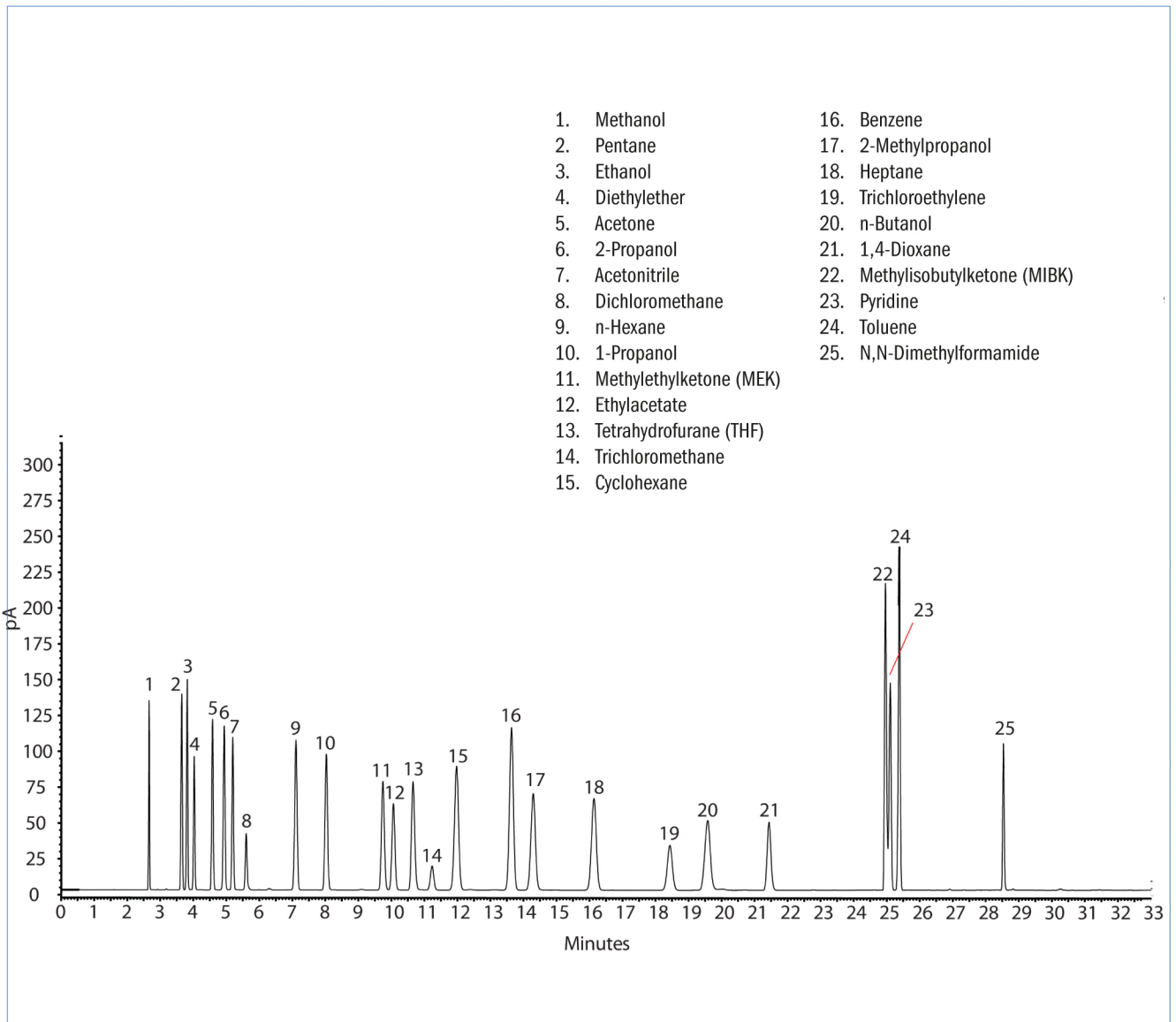

SAPIENS-624MS

www.teknokroma.es

Teknokroma Analitica S.A.- Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es