

**SAPIENS-624MS: EuPh Ethanol Impurities**

Column: SAPIENS-624MS: P/N, TR-591833


Dimensions: 30 m x 0,32 mm x 1,8 µm

Injection: 0.1 µl Test Alcohol Impurities (0.1%) 260°C, split: 1:20

Carrier gas: He, 10psi

Oven temperature: 40°C(12 min) to 240°C (5 min) @ 10°C/min

Detector: FID @ 260°C

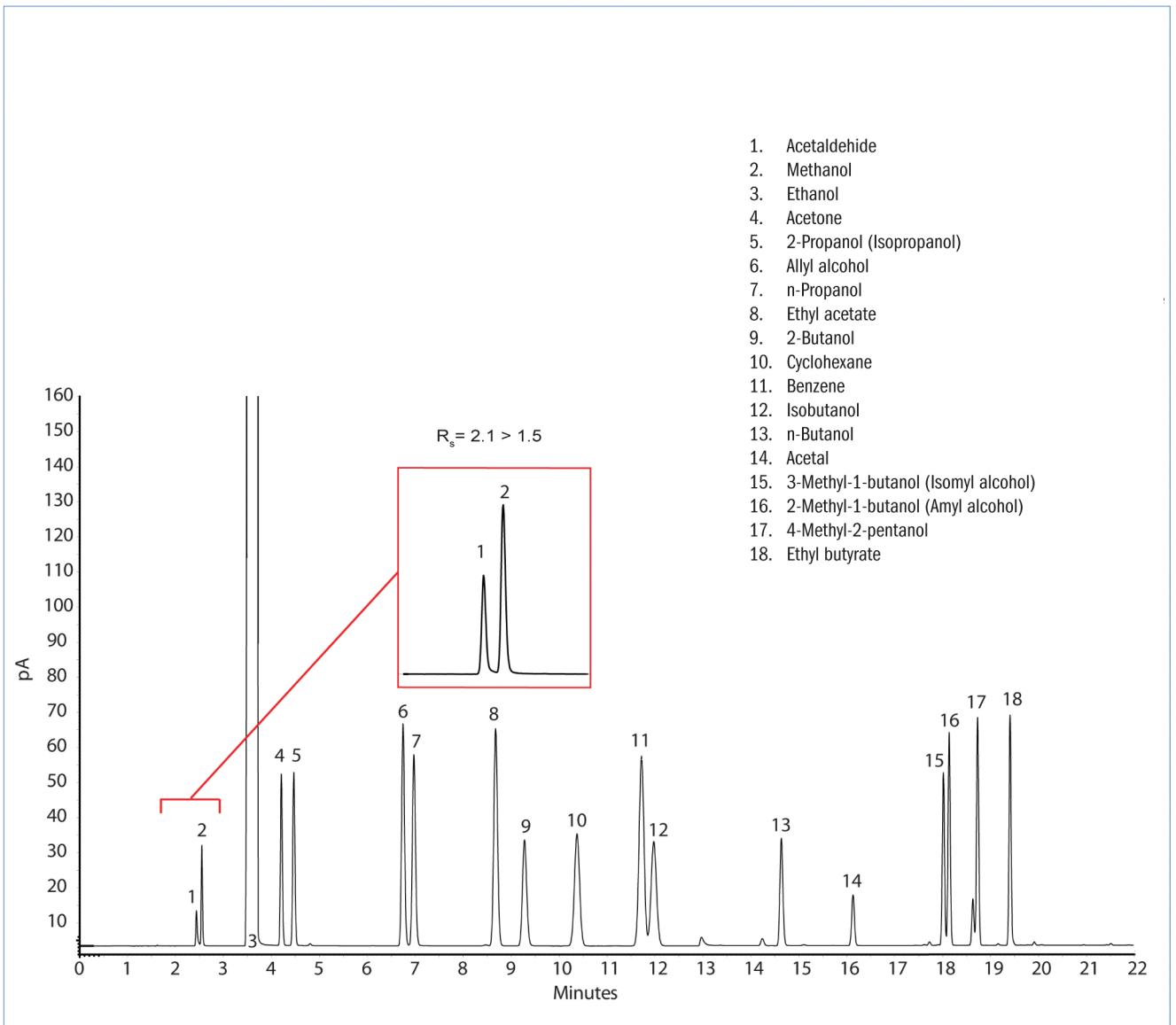

SAPIENS-624MS TKG 1305

[www.teknokroma.es](http://www.teknokroma.es)

Teknokroma Analitica S.A.- Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
 Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es