


## AMINES

Column: **TR-WAX.DB**, P/N TR-931035  
Dimensions: 30m x 0.53mm x 1.0µm  
Injection: 1 µL Amines mixture, Head Space, split 1:50, 260°C  
Carrier gas: H<sub>2</sub>, constant pressure 1.8 psi (12.40Kpa)  
Oven temperature: 60°C  
Detector: FID, 280°C

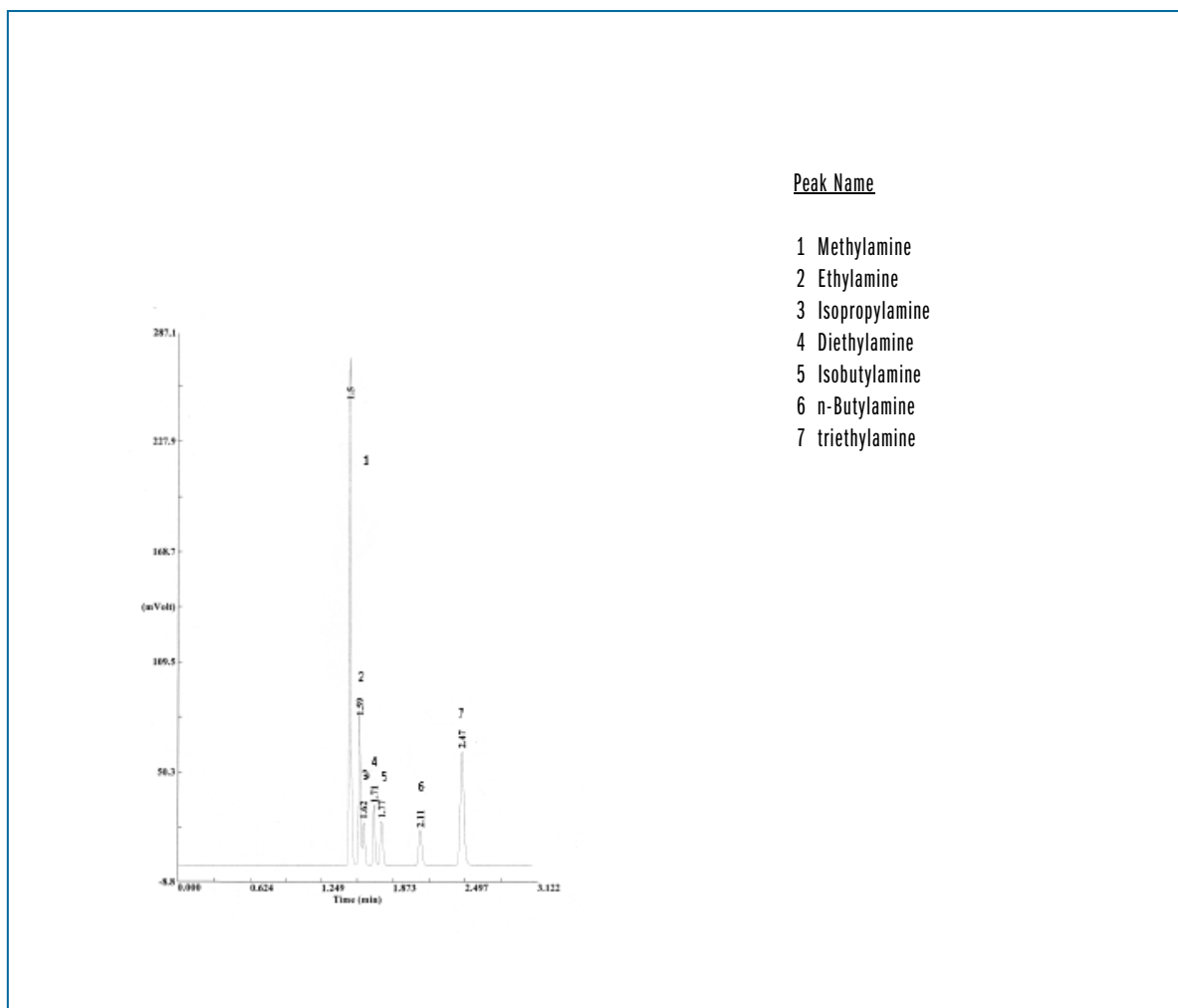

TKG 1058

[www.teknokroma.com](http://www.teknokroma.com)

Teknokroma S. Coop. C. Ltda. - Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain  
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es